




Protocol

Vena8 Endothelial+™ Biochip

Distributed in Europe by World Precision Instruments www.wpi-europe.com
UK Office +44 (0)1462 424700 wpiuk@wpi-europe.com
German Office +49 (0)6031 67708-0 wpide@wpi-europe.com

Vena8 Endothelial+ Biochip, Protocol #1: coating and cell seeding in Vena8 Endothelial+ biochips

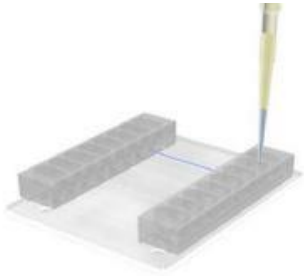

Step 1:

Cellix Vena8 Endothelial+ biochip is coated using a standard yellow tip pipette, by dispensing approximately 12 μL of protein (e.g. fibronectin, 100 $\mu\text{g}/\text{mL}$) into each microchannel. Note the excess of liquid on the entrance and exit ports.



Step 2:

The Vena8 Endothelial+ biochip is then placed in a humidified box which remains open for 1–1.5 hours in the CO_2 incubator. Alternatively, the biochip may be placed at 4°C for overnight coating.



Step 3:

After the incubation period, add 5 μL of $1.5 \times 10^6 / 100 \mu\text{L}$ ($\cong 15 \times 10^6 / \text{mL}$) of endothelial cells gently into each channel.

Note: concentration specified is for primary HUVEC.

The biochip is kept in the CO_2 incubator for 15–20 minutes for the cells to adhere. Observe the biochip under a microscope and top up all the reservoirs with 50 μL of media. Keep the biochip for 1.5–2 hrs in the CO_2 incubator.

Vena8 Endothelial+ Biochip Protocol #2: executing cell rolling, adhesion and migration assays under shear flow with Vena8 Endothelial+ biochips (manual version — not with VenaFlux platform)




Step 1:

Suspension cells (e.g. T cells, monocytes, platelets) are re-suspended in culture medium at an appropriate concentration (typically $2-5 \times 10^6/\text{mL}$) in an Eppendorf tube.




Step 2:

Using the Cellix Mirus Evo nanopump or the ExiGo pump, 10 μL of media is dispensed from pump output cable. Following this, the output cable is inserted into a specified channel on the Vena8 Endothelial+ biochip.

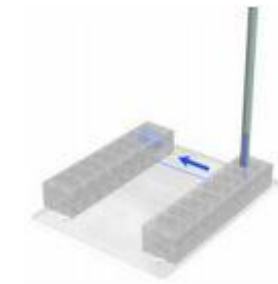


Step 3:

Then using the Cellix Mirus Evo nanopump or the ExiGo pump, 40 μL of the media is injected through the channel at a shear stress of 40 dynes/cm². This is done to wash the channel of cell debris. The waste is aspirated from the microwell of Vena8 Endothelial+ biochip with a pipette.

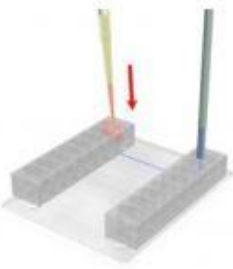


Step 4:

Cell sample is placed into the microwell of this channel on the Vena8 Endothelial+ biochip.

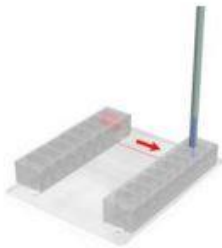

Step 5:

Cells are introduced into the channel, by specifying the desired shear stress on the FlowAssay software. The flow rate will be automatically calculated.

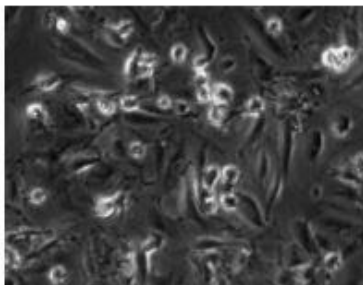

Step 6:

At each shear stress value, it is recommended that images of 3–5 fields of view of cell rolling and adhesion are acquired along the length of the channel.

Distributed in Europe by World Precision Instruments www.wpi-europe.com
UK Office +44 (0)1462 424700 wpiuk@wpi-europe.com
German Office +49 (0)6031 67708-0 wpide@wpi-europe.com