



WORLD
PRECISION
INSTRUMENTS

UV ROOM STERILIZER


CLEANRAY ROOM MODEL FOR DISINFECTION OF OFFICE, HOME OR LAB

DISINFECTION

99.99%
Effective



Finally, WPI offers an affordable way to disinfect your office, home or lab with the WPI CleanRay Room model. The CleanRay Room model stands upright with an aesthetic design to fit into any living space. A full 360° path to ensure sterilization of the entire room. To protect the occupants, the unit automatically shuts off when it detects movement in the room. It automatically restarts when the people are out of range. This prevents people and pets from coming into contact with the harmful UV rays.

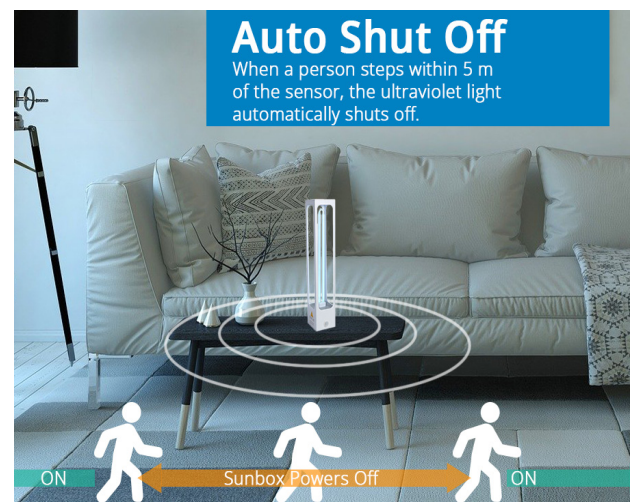


- **KILLS GERMS AND BACTERIA** - Hidden microorganisms lie waiting for you everywhere. With the WPI CleanRay Room model, UV light is all you need. Just set it up, turn it on and leave the room. Let the CleanRay remove your worries.
- **REDUCE YOUR CHANCE OF SICKNESS NOW!** Our hands are busy all day long. We touch hundreds of objects everyday, including our noses, eyes and mouths. The CleanRay disinfecting unit keeps you out of the doctor's office.
- **USE IT EVERYWHERE** - You will find yourself using this portable solution in every living space. Use it for kitchen counters, bathrooms, sitting areas, pet areas, offices, hotels and more.
- **100% SATISFACTION GUARANTEED** - Buy with confidence. If you are not satisfied, just send it back within 30 days for a full refund.
- **SANITIZES WITHOUT CHEMICALS OR IRRITANTS** - You can sanitize any surface (counters, bathroom fixtures, furniture, etc.) quickly and easily without chemical residue or harsh liquids or odors. Whether at home, work, dorm or bathroom, the CleanRay Room model is ideal for disinfecting almost everything you touch. It is not designed for use on people, pets or food.


UV ROOM STERILIZER

FEATURES:

- Combines two technologies in one unit - ozone and ultraviolet
- Quickly disinfect your whole room in 30 minutes (~430 ft²)
- Safety comes first with an automatic flashing light that indicates when the CleanRay is active
- Integral 360° degree induction motion sensor to shut off the light when someone enters the room
- Automatically restarts 10 seconds after everyone leaves the room
- Quick sterilization using 254 nm wavelength to breakdown RNA and DNA of microorganisms
- Stable base design
- Aesthetic design to place in any living space

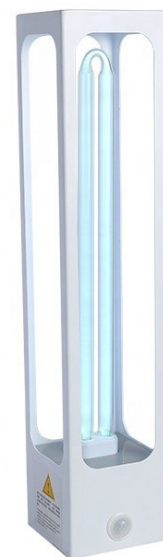

DUAL TECHNOLOGY USES OZONE AND ULTRAVIOLET

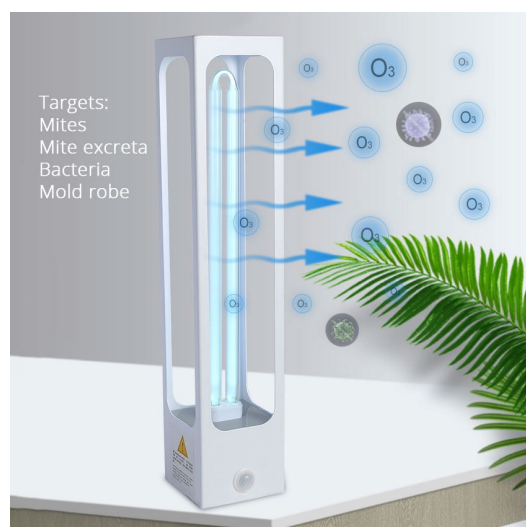


Ozone is a strong oxidant, which generates a large amount of free radicals and has a sterilizing and disinfecting effect.

The CleanRay Room model UV disinfection lamp is more effective than liquid disinfectant, without any irritating odor after disinfection. It is compact and easy to store, and it is safe and convenient to use.

The lamp bulb is made of a quartz material and emits 254 nm wavelength of light. The CleanRay has a long service life, high transmission rate and better sterilization effect.

SPECIFICATIONS

Power	38 W
Lamp Dimensions	100×100×482 mm (3.9×3.9×19")
Input Voltage	110 V/220 V
Material	Aluminum shell with Quartz lamp
Disinfection Method	UV and Ozone
Applicable Area	30-40 m ² (~430 ft ²)
Light Source	Quartz ultraviolet
Certificates	CE, RoHS
Control Mode	Remote control, time mode, motion sensor
Effective Disinfection	Bedroom: 15-30 min Kitchen/Bath: 15-30 min Living Room: 60 min

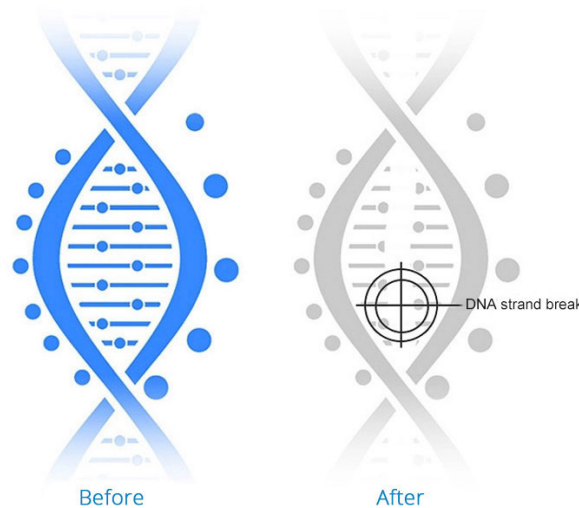

WORLD PRECISION INSTRUMENTS

USA: 175 Sarasota Center Boulevard, Sarasota FL 34240-9258 USA

Tel: (941) 371-1003 • **Fax:** (941) 377-5428 • **E-mail:** wpi@wpiinc.com • **Internet:** www.wpiinc.com

UK: 1 Hunting Gate, Hitchin, Hertfordshire SG4 0TJ England • Tel: 44 (0)1462 424700 • E-mail: wpiuk@wpi-europe.com

Germany: Pfingstweide 16, D-61169 Friedberg (Hessen) • Tel: +49 (0)6031 67708-0 • Fax: +49 (0)6031 67708-80 • E-mail: wpide@wpi-europe.com

China & Hong Kong: Rm 27A, No 8 Dongfang Rd., Pudong New District, Shanghai PRC • +86 6888 5517 • 400 688 5517 • ChinaSales@china.wpiinc.com

Brazil: Conselheiro Nabias, 756 sala2611, Santos-São Paulo 11045-002 Brazil • E-mail: info@brazil.wpiinc.com