


New Thermo Plate[®]

Microscope-Stage Automatic Thermocontrol System

TP series



Clear glass heating stage designed for specimen temperature control during the observation under microscopes.

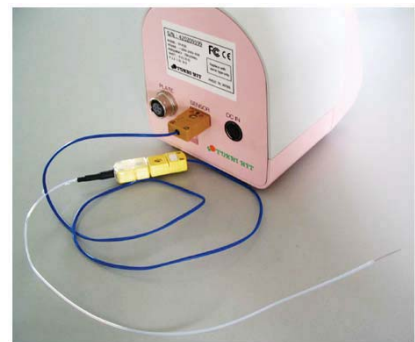

Stylish & Smart

Which is your color?



Feature 1 *Enables on site calibration*

A **sterilized sensor** measures the actual temperature. It allows you to calibrate yourself with easy operation and provide optimal environment for your laboratory. Furthermore, with the sensor, the system is able to be used as thermometer.





Heater temp. ⇄ Sensor temp.
【Monitor selector switch】



Sterilized sensor

Feature 2 *Stable Imaging, Effective for high magnification*

Continuous Current Control minimizes focus drift compared to the standard On-Off control method. Especially effective for observation when high magnification is required

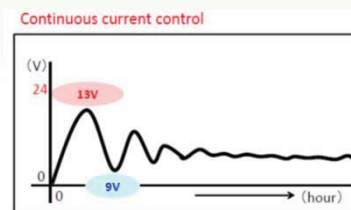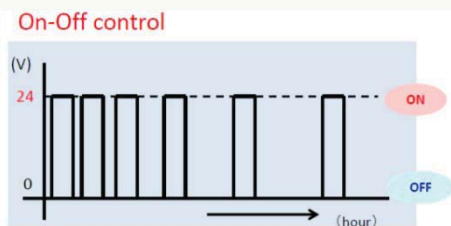

Feature 3 *Temperature management Software "TEM"*

Monitors, records, logs sample or plate temperature which the sensor measures. Save the data as CSV format which can be used on the spreadsheet.



Stylish & Smart
Which is your color ?

Size(Controller) :
105mm(W) x 116mm (H) x 148mm (D)



sales@wpi-europe.com